


08-815 Wireless High Power Headlight

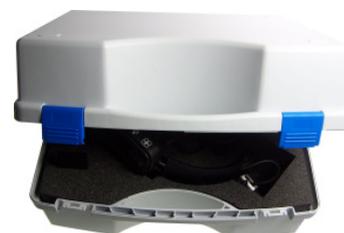

High performance, focusable headlight with white Power-LED. Wireless, light weight, comfortable headlight with two interchangeable Li-Ion battery packs at the headband. Especially suitable for long term use.

The 08-815 wireless LED-headlight provides an extremely **high illumination, long lasting operation** and a simple and **comfortable handling**. The **intensive illumination** is provided by a modern **high power LED**. The special high efficiency lenses focus the light to a **variable concentrated spot**. LED-headlights offer **the best possible price/performance relation** because of the typical long life and power efficiency of LED lamps and the modern Lithium-Ion technology. The batteries allow a continuous work of more than seven hours. The working time can be even extended by changing the batteries to fresh charged ones. Because of the **low total weight** and the ergonomic bolstered headband it is a **pleasure to wear** this headlight also for long times.

The light can be adjusted by a dimmer with "touch function". The ball joints make it possible to adjust the light beam exactly into any requested position for a restriction free optimal view.

The wall plug transformer for the charger includes connectors type EU, UK and USA.

Item code.	
08-815	Head light Set complete Delivery content: Headlight, 2x battery, Charger, wall plug power supply for charger, case
08-815-00	Head light (without batteries)
08-205-815	Battery pack
10-10815	Battery charger incl. power supply
08-811-40	Charging station without power supply
10-10212	wall plug power supply for charger (incl. plug adapter EU, US und UK)
08-815K	Storage case

Technical data:

Classification:	class I (annex IX, RL 93/42 MDD)
Classification:	1 M (EN 60825-1)
Battery:	Li-Ion
Power supply charger:	100-240V 50/60Hz
Charging time:	about 4 h if the battery is totally down
Working time:	> 7h (ca. 3,5 h/battery)
Working distance:	20...60 cm
Light spot:	focusable
Color temperature:	about 6.000 Kelvin
Weight:	326g with one battery; 376g with two batteries

Light intensity and Light spot size:

Working distance	Light spot Minimal Ø	Light spot Maximal Ø	Max. illumination [lux]
30cm	38mm	55mm	123200
40cm	53	85	71800
50cm	68	110	46800
60cm	80	135	31800