


## Integra<sup>®</sup> MicroFrance<sup>®</sup> Backbiting Forceps

Available in four sizes, the MicroFrance<sup>®</sup> Backbiting Forceps provide an outstanding cutting edge, thin shaft, and strong cutting capability while working in the nasal cavity.

Trying to use a shaver on the uncinate process can be difficult, so these specialized backbiters allow a reverse approach when removing or making an opening in the uncinate process during sinus surgery. They also provide a better alternative to traditional backbiting forceps because of their unique downward sloping shaft.

The slim profile and tapered shaft design fit ideally into the middle meatus to aid in reducing trauma to the middle turbinates during surgery. MCEN27-3 features a distinctive gentle upward curved shaft for accessing difficult-to-reach areas of the sinus.

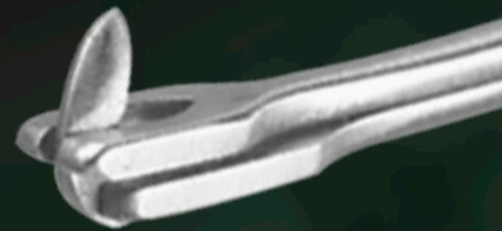

# Integra® MicroFrance®

## Backbiting Forceps

Reference	Description	Qty
MCEN27	Adult Forceps, 4.0 mm x 7.0 mm Bite	1
MCEN27-1	Adult Forceps, Fine, 3.5 mm x 7.0 mm Bite (not pictured)	1
MCEN27-2	Pediatric Forceps, 2.0 mm x 4.0 mm Bite (not pictured)	1
MCEN27-3	Pediatric Mini-Forceps, 2.0 mm x 3.5 mm Bite	1


MCEN27-3



MCEN27



MCEN27-3

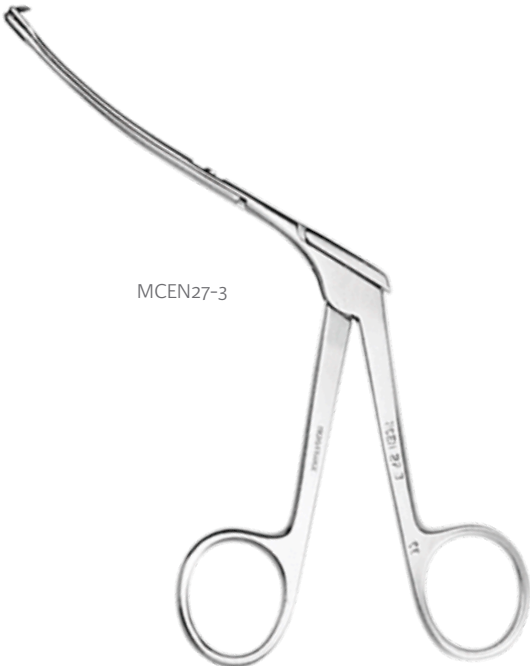

MCEN27



The MicroFrance® brand is a premier line of ENT instruments, designed with world-renowned ENT surgeons to facilitate new techniques and improve patient outcomes.

Availability of these products might vary from a given country or region to another, as a result of specific local regulatory approval or clearance requirements for sale in such country or region.

- Always refer to the appropriate instructions for use for complete clinical instructions.
- Non contractual document. The manufacturer reserves the right, without prior notice, to modify the products in order to improve their quality.
- Warning: Applicable laws restrict these products to sale by or on the order of a physician.

For more information or to place an order, please contact:

United States, Canada, Asia, Pacific, Latin America

USA 800-654-2873 ▪ 888-980-7742 fax

International +1 609-936-5400 ▪ +1 609-750-4259 fax

[integralife.com/contact](http://integralife.com/contact)

Manufacturer:



Integra MicroFrance SAS  
Le Pavillon, 03160  
Saint-Aubin-le-Monial ▪ France