


Integra[®] MicroFrance[®]

Duncavage Sphenoid Sinus Instrument Set

Developed in conjunction with James A. Duncavage, MD

The Duncavage Sphenoid Sinus Set was developed to allow access to the interior of the sphenoid sinus. These instruments help surgeons remove cysts and bone, as well as strip mucosa and palpate masses within the sphenoid sinus.

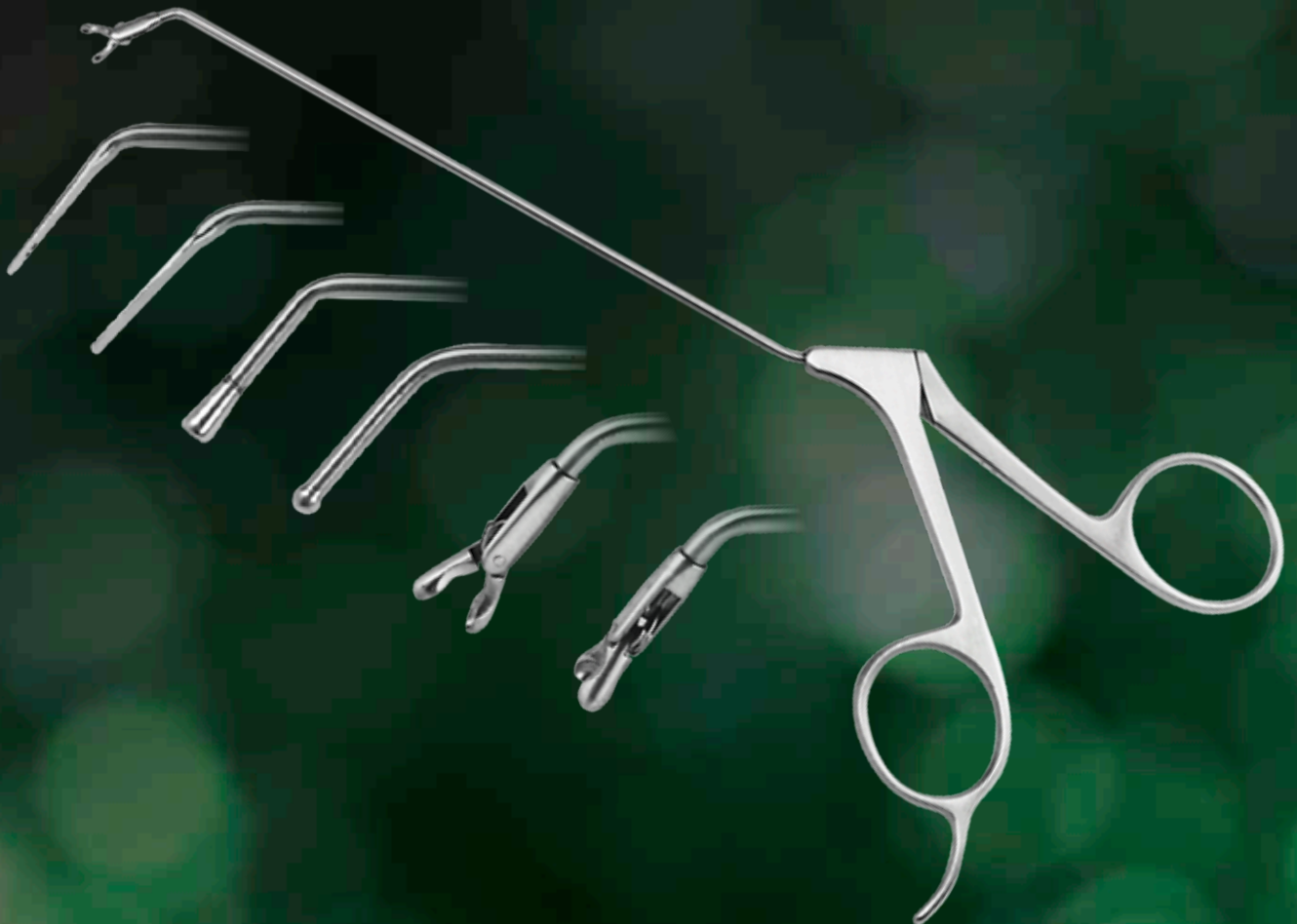

Integra® MicroFrance®

Duncavage Sphenoid Sinus Instrument Set

All of the Duncavage instruments feature a 25° bend superiorly and a 25° bend inferiorly at the distal tips that allow the surgeon to help access the sphenoid sinus. They also feature reverse-action handles that help reduce the amount of movement at the distal tips when actuated.

Ordering Information

Product	Description	Quantity
MCEN150*	Duncavage Sphenoid Sinus Set	1 set
Set Includes:		
MCEN151	Sphenoid Sinus Vertical Cup Forceps, 2.5mm, 94mm (3¾")	1
MCEN152	Sphenoid Sinus Horizontal Cup Forceps, 7mm, 93mm (3¾")	1
MCEN153	Sphenoid Sinus Vertical Thru-Cut Forceps, 3mm, 93mm (3¾")	1
MCEN154	Sphenoid Sinus Horizontal Thru-Cup Forceps, 7.3mm, 93mm (3¾")	1
MCEN155	Sphenoid Sinus Blunt Probe, 2.5mm, 140mm (5½")	1
MCEN156	Sphenoid Sinus Blunt Probe w/Suction, 2.9mm, 131mm (5¼")	1
MCEN157	Sphenoid Sinus Freer Elevator, 3 x 17mm, 140mm (5½")	1
MCEN158	Sphenoid Sinus Freer Elevator w/Suction, 2.7 x 12mm, 132mm (5¼")	1
3717002	Instrument Tray, 52cm L x 25cm W x 2.8cm D, (20½" L x 9¾" W x 1⅛" D)	1

*This set is only available in the US. If you are not in the US, please order the above items individually.



The MicroFrance® brand is a premier line of ENT instruments, designed with world-renowned ENT surgeons to facilitate new techniques and improve patient outcomes.

Availability of these products might vary from a given country or region to another, as a result of specific local regulatory approval or clearance requirements for sale in such country or region.

- Non contractual document. The manufacturer reserves the right, without prior notice, to modify the products in order to improve their quality.
- Warning: Applicable laws restrict these products to sale by or on the order of a physician.
- Consult product labels and inserts for any indication, contraindications, hazards, warnings, precautions, and instructions for use.

For more information or to place an order, please contact:

United States, Canada, Asia, Pacific, Latin America

USA 800-654-2873 ▪ 888-980-7742 fax

International +1 609-936-5400 ▪ +1 609-750-4259 fax

integralife.com/contact



Integra MicroFrance SAS

Le Pavillon, 03160

Saint-Aubin-le-Monial ▪ France