

## Integra<sup>®</sup> MicroFrance<sup>®</sup>

### Levine Orbital Decompression Instrument Set

Developed in conjunction with Howard L. Levine, MD

The Levine Orbital Decompression Instrument Set is suited for intranasal endoscopic orbital decompression. The Levine Perforator is designed to create the initial opening through the medial orbital wall (lamina papyracea). The 15 degree angled suction curettes and picks are designed to help gently remove the lamina. The angled sickle knife incises the periorbital to liberate the orbital fat. The serrated jaws of the angled forceps do not completely close, allowing for a firm, yet atraumatic grip to assist in teasing the orbital fat into the sinus cavity. The 10-instrument set is provided with an instrument tray.





# Integra® MicroFrance®

## Levine Orbital Decompression Instrument Set

### Ordering Information

Product	Description	Quantity
MCEN160*	Levine Orbital Decompression Instrument Set	1 set
<b>Set Includes:</b>		
MCEN161	Orbital Ring with Suction, 208mm (8¼")	1
MCEN162	Orbital Knife, 185mm (7¼")	1
MCEN163	Orbital Suction Pick, 90°, 208mm (8¼")	1
MCEN164	Orbital Suction Pick, 45°, 208mm (8¼")	1
MCEN165	Orbital Suction Pick, 135°, 208mm (8¼")	1
MCEN167	Orbital Suction Curette, 90°, 209mm (8¼")	1
MCEN168	Orbital Sickle Knife, 230mm (9")	1
MCEN169	Orbital Forceps, 45°, 118mm (4⅝")	1
MCEN170	Orbital Curette Dissector, 250mm (9¾")	1
MCEN171	Orbital Perforator Elevator, 250mm (9¾")	1
3717002	Instrument Tray, 52cm L x 25cm W x 2.8cm D, (20½" L x 9¾" W x 1⅛" D)	1

\*This set is only available in the US.  
If you are not in the US, please order the above items individually.

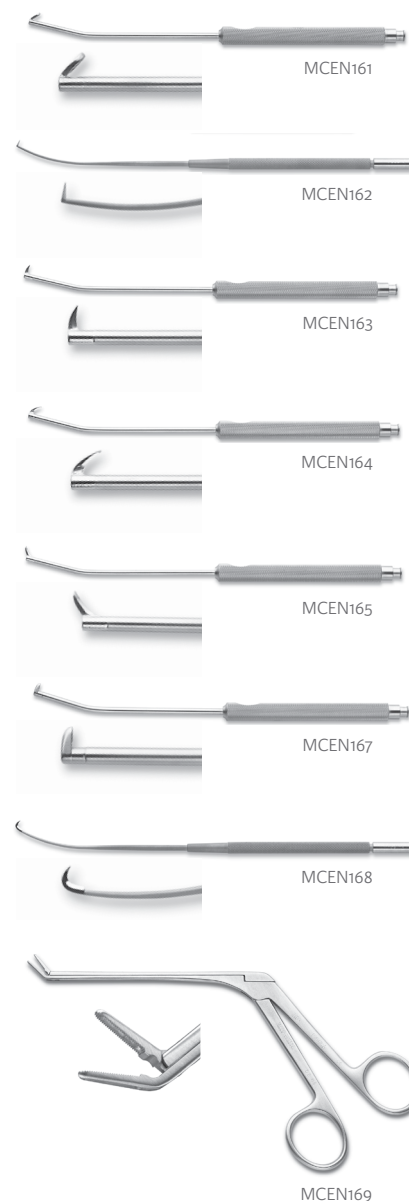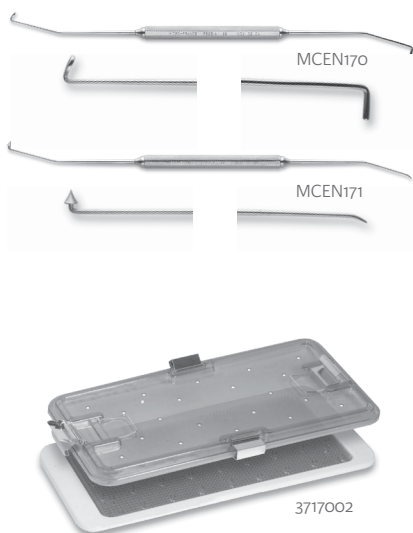

The MicroFrance® brand is a premier line of ENT instruments, designed with world-renowned ENT surgeons to facilitate new techniques and improve patient outcomes.

- Availability of these products might vary from a given country or region to another, as a result of specific local regulatory approval or clearance requirements for sale in such country or region.
- Non contractual document. The manufacturer reserves the right, without prior notice, to modify the products in order to improve their quality.
- Warning: Applicable laws restrict these products to sale by or on the order of a physician.
- Consult product labels and inserts for any indication, contraindications, hazards, warnings, precautions, and instructions for use.

For more information or to place an order, please contact:

United States, Canada, Asia, Pacific, Latin America

USA 800-654-2873 ▪ 888-980-7742 fax

International +1 609-936-5400 ▪ +1 609-750-4259 fax

[integralife.com/contact](http://integralife.com/contact)

