

Integra®

Ruggles®-Redmond™ Distraction Screws


INTEGRA®
LIMIT UNCERTAINTY



Limit Uncertainty with symmetrical distraction of vertebral bodies for unobstructed views of disc space.

Used with Anterior Cervical Fusion (ACF) Systems

- Packaged individually in 12, 14, and 16mm sizes
- Available sterile and non-sterile
- For single use only
- Sizes clearly marked on each distraction screw
- Contains a sharp self-tapping tip design for insertion
- Compatible with Caspar Type Screwdrivers

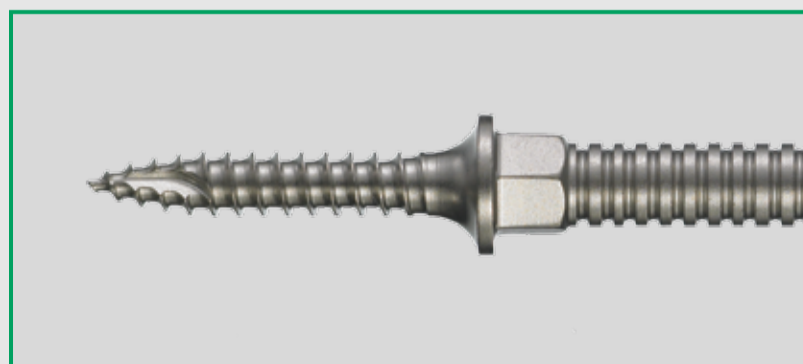


- The Ruggles-Redmond **ACF** Instrumentation System facilitates exposure of the operative site with minimal dissection of and damage to muscle tissue.
- The Ruggles-Redmond Distraction system permits symmetrical distraction of vertebral bodies for unobstructed view of disc space and easier graft placement.
- Distraction from the center of the vertebral body reduces limitations caused by narrow disc space
- For single level or multi level **ACF** procedures
- Available in left and right directions

Distraction Screws are used in conjunction with orthopedic retractors for distraction of vertebral bodies during cervical procedures.

Distraction Screws

Reference	Description
R6367A	12mm Sterile Distraction Screw
R6398A	14mm Sterile Distraction Screw
R6397A	16mm Sterile Distraction Screw
RN2877	12mm Non-Sterile Distraction Screw
RN2878	14mm Non-Sterile Distraction Screw
RN2879	16mm Non- Sterile Distraction Screw

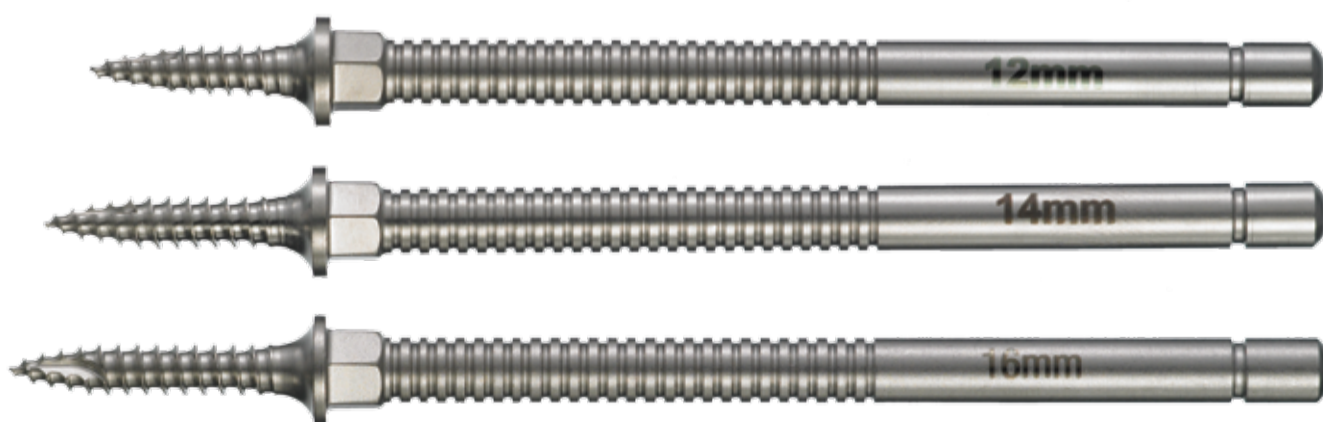


Clinical Applications

- Ruggles-Redmond Anterior Cervical Fusion (ACF) instrumentation is used to treat disc disease, pseudarthrosis, cervical spine injury and tumors involving vertebral bodies.

ACF Instrumentation System

Reference	Description
1. R6318-R6319	ACF Bone Grafting Forceps
2. R6317	Caliper
3. R6314-R6316	Bone Graft Impactor
4. R6396	Caspar Type Screwdriver
5. R9395	Caspar Type Twist Drill
6. R6388 or R6393	Caspar Type Drill Guide
7. RB6387 – RB6394	Caspar Type Distractor
8. R6312	ACF Transverse Retractor Body
9. R6555 – R559	ACF Medial Retractor Blades
10. R6560 – R6564	ACF Lateral Retractor Blades
11. R6310	ACF Blade Ejector
12. R6313	ACF Longitudinal Retractor Body
13. R6320 or R6321	ACF Fenestrated Illiac Blades
14. R6550 or R6554	ACF Longitudinal Retractor Blades

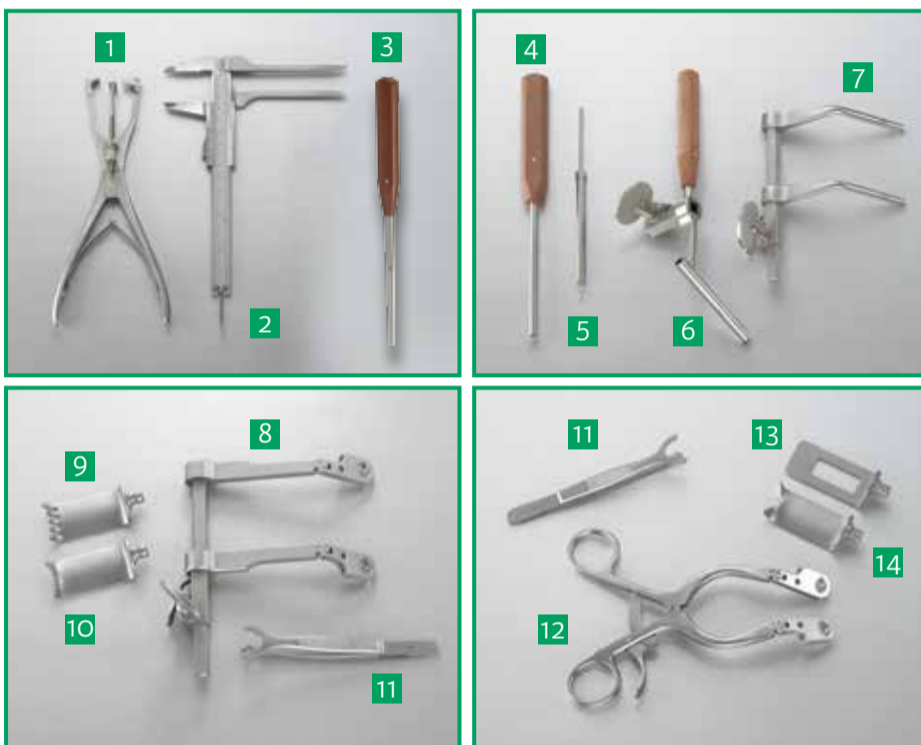

For more information or to place an order, please contact:

Integra ■ 589 Davies Drive, York, PA 17402

877-468-5572 USA ■ 914-789-7094 outside USA

609-750-4259 fax ■ 914-592-8056 outside USA fax

integralife.com

Ruggles, Integra and the Integra logo are registered trademarks of Integra LifeSciences Corporation in the United States and/or other countries. Redmond is a trademark of Integra LifeSciences Corporation.

©2012 Integra LifeSciences Corporation. All rights reserved. PM12113